



Guided Start Drill
37152*



Guided Twist Step Drill
Ø2.0 Pilot (#2) 32746**
Ø2.4/2.8 (#3) 35844**
Ø2.8/3.2 (#4) 35875**
Ø3.2/3.6 (#5) 35876**
Ø3.8/4.2 (#6) 35877**
Ø4.2/4.6 (#7) 38045**
Ø4.2/5.0 (#8) 37945**



Prosthetic Adapter
29167



Surgical Torque
Wrench Manual
28839



Guided Twist Drill
Ø1.5 x 20mm
33066*



Guided Anchor
Pin Ø1.5 mm
30909



Guided Anchor Pin
Ø1.5 mm short shaft
34761*



Guided Tissue Punch
NP 38059
RP 37154
WP/6.0 (5.0) 37155*
WP/6.0 (5.5) 38060*



Nobel
Biocare™

NobelParallel™ CC
Guided PureSet™
301163

NobelParallel™ CC Guided PureSet™

Anchor Pin 800 rpm max

Tissue Punch 800 rpm max

Drill Guide 800 rpm max

Cortical Drill 800 rpm max

Screw Tap 45 Ncm max 25 rpm

Implant Mount 45 Ncm max 25 rpm

Template Abutment

NP Ø3.75

RP Ø4.3

RP Ø5.0

WP Ø5.5

P Guided

Bone density Soft Medium Dense

NP Ø3.75 2 [3] 2 3 C [S] 2 3 4 C S

RP Ø4.3 2 3 [5] 2 3 5 C [S] 2 3 5 C S

RP Ø5.0 2 3 5 [6] 2 3 5 6 C [S] 2 3 5 6 C S

WP Ø5.5 2 3 5 7 [8] 2 3 5 8 C [S] 2 3 5 8 C S

[-] = optional

C = cortical drill

S = screw tap

Handpiece Connector

Screwdriver Unigrip




Guided Drill Guide NP
#2: Article 32814
#3: Article 32817
#4: Article 35882



Guided Drill Guide RP 4.3
#3: Article 32818
#5: Article 35883



Guided Drill Guide RP 5.0
#2: Article 32816*
#3: Article 38046*
#5: Article 35885*
#6: Article 32825*

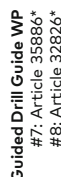

Guided Drill Guide WP
#7: Article 35886*
#8: Article 32826*



Guided Cortical Drill
NP Ø3.75 38061*
RP Ø4.3 38062*
RP Ø5.0 38063*
WP Ø5.5 38064*



Guided Screw Tap
7-10 mm
NP Ø3.75 38047*
RP Ø4.3 38050*
RP Ø5.0 38053*
WP Ø5.5 38056*



Guided Screw Tap
11.5-13 mm
NP Ø3.75 38048*
RP Ø4.3 38051*
RP Ø5.0 38054*
WP Ø5.5 38057*




Guided Screw Tap
15-18 mm
NP Ø3.75 38049*
RP Ø4.3 38052*
RP Ø5.0 38055*
11.5-15 mm
WP Ø5.5 38058*

* Article not included in this kit

** Article not included in this kit and also available in other lengths

*** Screwdrivers available in other lengths

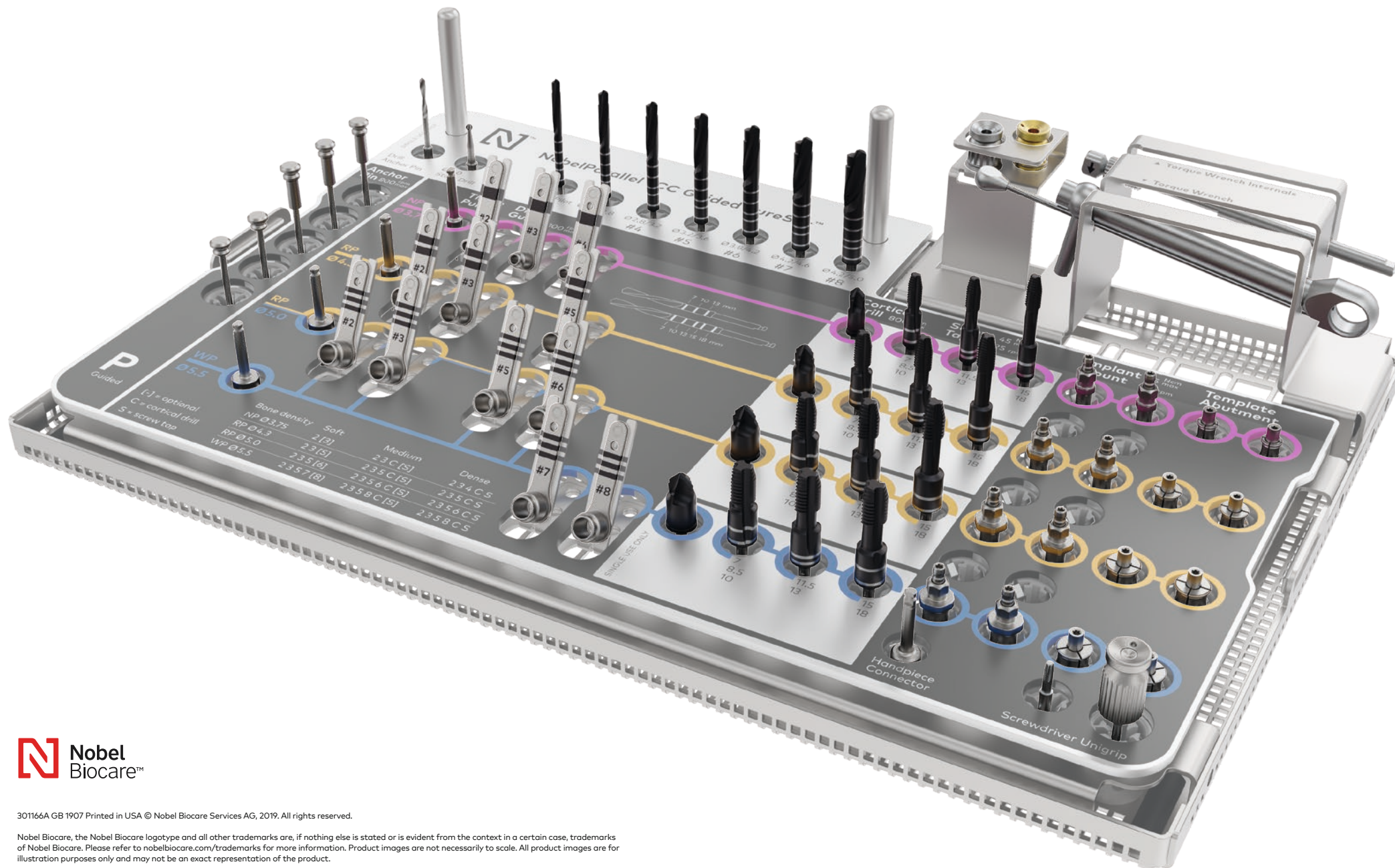

301166A GB 1907 Printed in USA © Nobel Biocare Services AG, 2019. All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Please refer to [nobelbiocare.com/trademarks](https://www.nobelbiocare.com/trademarks) for more information. Product images are not necessarily to scale. All product images are for illustration purposes only and may not be an exact representation of the product.

Disclaimer: Some products may not be regulatory cleared/released for sale in all markets. Please contact the local Nobel Biocare sales office for current product assortment and availability. For prescription use only. Caution: Federal (United States) law restricts this device to sale by or on the order of a clinician, medical professional or physician. See Instructions for Use for full prescribing information, including indications, contraindications, warnings and precautions.